


# The World of Dental Art

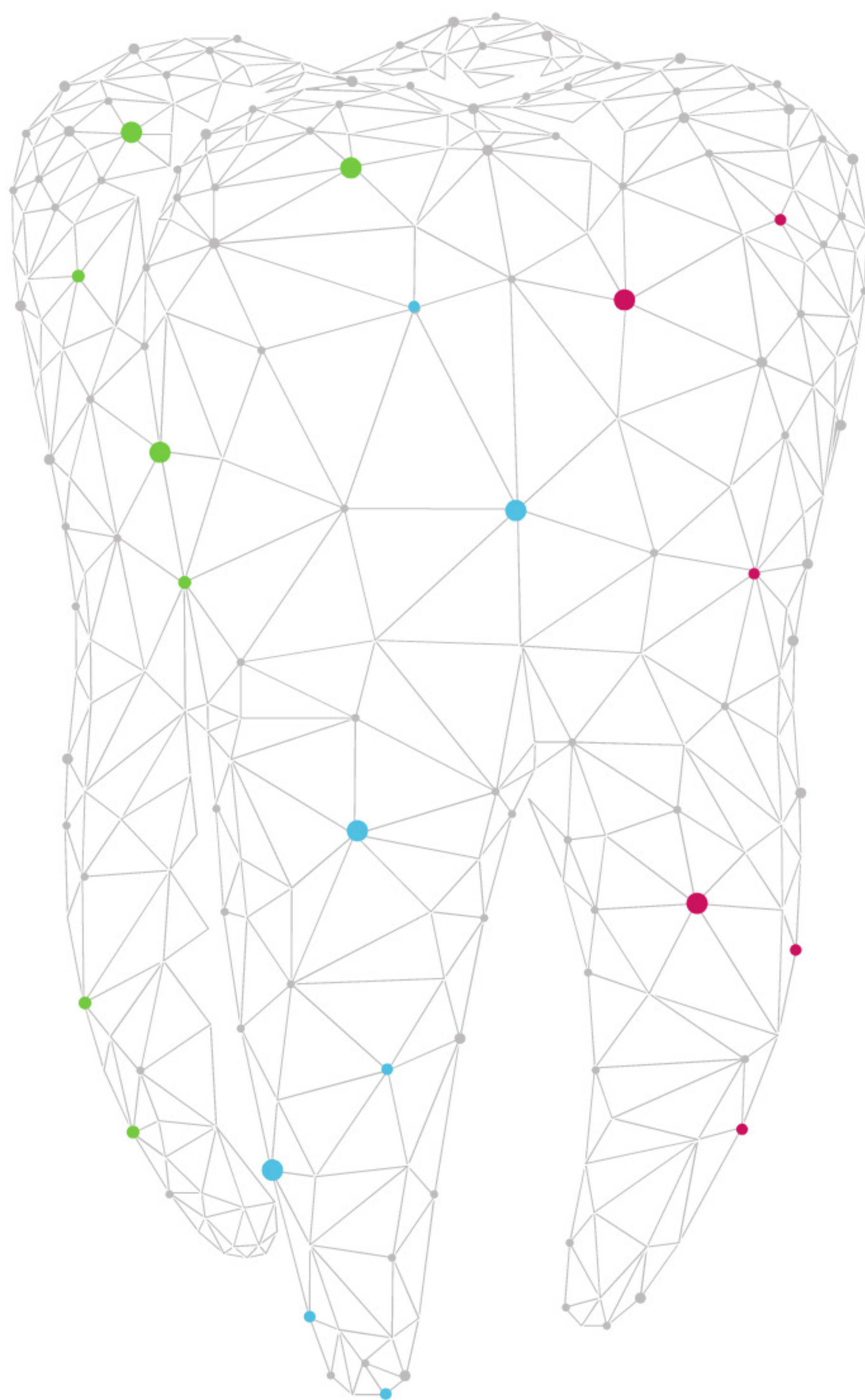

**PRODUCT  
CATALOGUE**

## DS Technology LTD

### INNOVATION THAT MAKES THE DIFFERENCE

DS Technology is a leading research and development company for Dental products with innovations, excellent quality and service to simplify the working process in dental lab and clinics.

DS Technology developing a range of dental solutions requires talent, experience, and an eye for detail. Headed by a skilled dental technician with 24 years of hands-on experience.

We are proud to provide efficient, professional and aesthetic solutions for dental labs, dentistry professional experts and their assisting teams globally. Our specialized products and technologies have made procedures more efficient, accurate and practical for dental professionals, as well as more comfortable and excellent aesthetic restorations with better fit and finish for their patients

We offer a self-cure revolutionary new material developed specifically to stain and glaze for composite, acrylic and PMMA type materials designed especially for precision customization and natural look, beautiful and esthetic restorations. Each material improves work efficiency, enhances results and increases endurance, resulting in a highly reliable application.

Our original, advanced thermoplastic denture base materials support the manufacturing of high standard for aesthetic removable dentures in today's modern dentistry.

DS Technology will continuously strive to supply our customers with innovative and improved products while simultaneously contributing to the field of dentistry and providing our customers with products of outstanding quality.

DS Technology LTD quality standards also meet the stringent specifications of the ISO 13485 and the CE mark.



# VITOSIL®

VPS Impression system  
for any technique

**Precision at it's best!**



- High hydrophilicity
- Highest precision
- Excellent flow properties
- Excellent reproduction of details
- Very low shrinkage

- Long term dimensional stability
- Optimal final hardness, Tear resistant
- Long working time
- Short intra oral setting time
- Simplicity and accuracy



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## Table of Contents

<b>Denture</b>	<b>4+5</b>	<b>Die Preparation</b>	<b>32</b>
Top Natur Glaze	6	Top PC Spacer	33
VITO CRYL Seal - IQ®	6	Top PC Spacer - PCS ISO	33
Top Airstain Glaze	7	Top PC Spacer - KITS	34
Top Airstain Glaze - KITS	8	Top Die Spacer	35
Top Air Opaque	9	Top Die Lube	35
Top Paqubond	9	Top Die Lube - Pen Applicator	35
Top Air Opaque - KITS	10	Top Rubbsep	36
		Top Blockout Wax	36
<b>DSTVITO – FLEXI DENTURE MATERIALS</b>	<b>11+12</b>	<b>CAD/CAM</b>	<b>37+38</b>
DSTVITOflexi	13	Top Scan Wax	39
DSTVITOacry	14	Bari x-ray	39
DSTVITOlux-f	15		
DSTVITOkeep	16	<b>Wax Product</b>	<b>40</b>
DSTVITOpplast	17+18	<b>Dipping Waxes</b>	
DSTVITOcetal	19	Top Dip Wax S	41
DST CLEAN & RENEW	20	Top Dip Wax M	41
DST LUSTERLIQ	20	Top Dip Wax H	41
DST LUBRILIQ	20	<b>Crown and Bridge Waxes</b>	
DST ISOLIQ F	21	Top Modeling Wax	42
DST CEM F	21	Top PC Modeling Wax	42
DST FUSE FLEX	21	Top Sculpturing Wax	42
Top Isoliq A	22	<b>Special Waxes</b>	
Top Gyps Rid	22	Top Diagnostic wax	43
Top Bite Rims	22	Top Cervical Wax	43
		Top Scan wax	44
<b>Ceramic Products</b>	<b>23+24</b>	Top Blockout Wax	44
Top Ceraliq	25	Top Wax Finish	44
Top Cera Glazer	25	Uni-Wax Adhesive	44
Top Marker	25	<b>Sticky Waxes</b>	
Top Cerafix Gel	26	Top Sticky Wax	45
Top Cera Contour	26	Top Spec Sticky Wax	45
Top Cera Sep	27		
Top Cera Sep - Pen Applicator	27		
Top Cera Sep Spray	27		
<b>Model Production</b>	<b>28</b>		
TOPDEB Debubbler	29		
Top Vacudeb	29		
Top Sep	29		
Top Smart Sep	30		
Top shield&Shine Spray	30		
Top Blockout Wax	30		
Top Die Hardener& Sealer	31		
Top Die Hardener& Sealer - Pen Applicator	31		
Top Die Sealer	31		



# Denture Products

Top Natur Glaze | VITO CRYL Seal - IQ® | Top Airstain Glaze  
Top Airstain Glaze - KITS | Top Air Opaque  
Top Paqubond | Top Air Opaque - KITS



Dental Materials Production  
[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## Top Natur Glaze

Clear acrylic liquid for surface sealing of a denture and any acrylic appliance .

Brilliant natural finish, without mechanical polishing.

Top Natur Glaze produces a natural wet-looking protective glaze over permanent and provisional acrylic allowing the base color to show through. Self-curing, increases color stability. Fills micro voids to reduce bacteria and odor problems.

Bonds to all dental resins, including methyl and ethyl methacrylates, polycarbonates, bis-acryl resins, resin crowns and laminates, vinyl methacrylate copolymers, acrylic denture teeth denture base, and composite resins.



- | Brilliant natural finish, without mechanical polishing.
- | Fully transparent, no additional polishing required.
- | Creates a very smooth surface.
- | High gloss, colourless, strong bonding
- | Also protects against plaque deposits and discolouration.
- | Applicable for composite restoration.
- | Rapid air dry within 20 seconds.

Part num	Color	Amount
26-272	Transparent	1/2 oz
Thinner		
26-275	Transparent	1/2 oz

## VITO-CRYL Seal - IQ®

Clear acrylic spray to seal & shine the tissue surface a denture and any acrylic appliance

Clear acrylic spray to seal&shine a fills open pores the tissue side of a denture and any acrylic appliance, to seal the tissue surfaces from both bacteria and plaque and Brilliant natural finish, without rough spots.



- | EASY TO USE!
- | Impressive results in less time
- | Micro thin layer
- | Creates an excellent Seal&Shine smooth surface
- | Improve patient comfort,
- | Dries quickly - 60 to 90 seconds
- | Chemically bonds to all-types of acrylic resin

Part num	Color	Amount
18-101	Transparent	1 oz
18-102	Transparent	2 oz



Top Airstain Glaze



**TOP AIRSTAIN GLAZE**  
Self-Cure Colored Acrylic Liquid Resin System

**PRECISION CUSTOMIZATION**  
Natural look restorations

**For aesthetic color modification of any acrylic, composite and PMMA type materials.**

Top Air Stain Glaze colors is a next generation self-cure acrylic liquid stain a new, specially developed to stain and glaze composite, acrylic and PMMA type materials.

Top Air Stain Glaze is an self-cure acrylic stain, creating a natural effects and glaze finish for quick characterizing, shading and staining of provisional and permanent acrylic and composite restorations - artificial teeth, acrylic dentures, implant retained dentures, composite crowns and bridge, temporary resin crowns and bridges and acrylic permanently fixed with implants within seconds.

Top Air Stain Glaze are available in an 11 natural colors for perfect blending and matching. All 11 Art stains can be applied alone or blended together to achieve different effects, available for color characterization and the shade effects of natural teeth and the pigmentation of the gums.

Translucent and stable, these stains do not change tint or hue during application or over time. It dry quickly in less than 20 seconds, and the color can be built up gradually.

- | Beautiful and esthetic restoration
- | Immediate self-curing for results within seconds
- | Powerful properties bond and cure to all dental resins
- | Applicable for acrylic, composite and PMMA type materials
- | Creating effects and glaze in all resin type prosthetics
- | Personalized, precision color matching
- | Translucent and stable
- | Effects of natural realistic appearance

<b>Part no</b>	267000	262701	262702	262703	262704	262705	262706	262707
<b>Color</b>	Violet	White	Brown	Gray	Blue	Pink	Dark Pink	Yellow
<b>Amount</b>	10ml	10ml	10ml	10ml	10ml	10ml	10ml	10ml
<b>Part no</b>	262708	262709	262710	26275	26274			
<b>Color</b>	Vanilla A+B	Beige C+D	Orange	Thinner	Brush Cleaning			
<b>Amount</b>	10ml	10ml	10ml	1/2 oz	1/2 oz			



Top Airstain Glaze - KITS



**Professional Kit - 9 Colors**

**Part num - 26- 270**

**9 colors in 10ml Bottle - Pink, Dark Pink, Yellow, Brown, Gray, Blue, White and two cervical colors: Vanilla- A+B, Beige- C+D**

**Top Natur Glaze - clear - 1/2 oz**

**Thinner - 1/2 oz**

**Brush Cleaning Thinner - 1/2 oz**

**1 Mixing Well**

**2 Brush**

**Tooth Aesthetic Kit - 4 Colors**

**Part num - 26- 273**

**4 colors in 10ml Bottle - Gray, Yellow, Vanilla- A+B, Beige- C+D**

**Top Natur Glaze - clear - 1/2 oz**

**Thinner - 1/2 oz**

**Brush Cleaning Thinner - 1/2 oz**

**1 Mixing Well**

**2 Brush**

**Denture Aesthetic Kit - 4 Colors**

**Part num - 26- 271**

**4 colors in 10ml Bottle - Pink, Dark Pink, Blue, White**

**Top Natur Glaze - clear - 1/2 oz**

**Thinner - 1/2 oz**

**Brush Cleaning Thinner - 1/2 oz**

**1 Mixing Well**

**2 Brush**



**Dental Materials Production**

www.dst-dental.com | Email: info@dst-dental.com

Top Air Opaque



**Top Air Opaque**  
Self-cure acrylic opaque liquid - metal mask

**Now it is an easy task**  
**Successfully block the dark metal hue**

Top Air Opaque - is a ready to use, self-cure metal mask acrylic opaque liquid, perfectly covers the metal surface in a very thin coat. Top Air Opaque is a quick drying, compact colored acrylic opaque system of dentin & gingiva shades for use in dental restorations where it is needed to mask metal as denture base framework, veneers, copings, bridgework and implant bar for pleasing aesthetics results. A rich variety of colors teeth and gums, allow an easy and comfortable fit for all restorations.

**Note: Not suitable for mask areas that are not covered by acrylic - such as clasps.**

- | Ready to use, Self-cure
- | Easy and clean application
- | Super-thin coating
- | Fast 30 seconds-drying time
- | Compatible with all acrylic & composite
- | Wide range of colors
- | Aesthetic restorations
- | Affordable, cost effective

Part no	262901	262902	262903	262904	262905	262906	26295
Color	Pink	Light Pink	Dark Pink	Snow	Craem A+B	Beige C+D	Thinner
Amount	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz

Top Paqubond



**TOP PAQUBOND P/L**

Top Paqubond is a complement product of Top Air Opaque - is intended to apply as a first coat before applying Top Air Opaque and his purpose a increasing adhesion.

The system comprises Powder and Liquid.

Top Paqubond is specially development to strengthen the adhesion of the Top Air Opaque. Therefore the use of two products - Top Paqubond with Top Air Opaque together - according to the instructions for use - this to enhance bonding properties.

**Note: Not suitable for mask areas that are not covered by acrylic - such as clasps.**

	<b>POWDER</b>	<b>LIQUID</b>	<b>SET</b>
Part no	26-2922	26-2921	26-2923
Amount	7g	10ml	

Top Air Opaque - KITS



**COMPLETE KIT**

**Part num - 26-290**  
**6 colors in 1/2 oz Bottle - Pink, Light Pink, Dark Pink, Snow, Cream - A+B, Beige - C+D**  
**Top Air Opaque Thinner - Transparent - 1/2 oz**  
**Top Paqubond- Special to Top Air Opaque - POWDER 7g**  
**Top Paqubond- Special to Top Air Opaque - LIQUID 10ml**  
**Mixing Pad**  
**Micro Applicator Brushes**

**STARTER KIT**

**Part num - 26-291**  
**3 colors in 1/2 oz Bottle - Pink, Snow, Cream - A+B,**  
**Top Air Opaque Thinner - Transparent - 1/2 oz**  
**Top Paqubond- Special to Top Air Opaque - POWDER 7g**  
**Top Paqubond- Special to Top Air Opaque - LIQUID 10ml**  
**Mixing Pad**  
**Micro Applicator Brushes**





Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)



# DSTVITO – FLEXI DENTURE MATERIALS

DSTVITOflexi | DSTVITOacry | DSTVITOlux-f | DSTVITOkeep  
DSTVITOpplast | DSTVITOCetal | DST CLEAN & RENEW | DST LUSTERLIQ  
DST LUBRILIQ | DST ISOLIQ F | DST CEM F | DST FUSE FLEX  
Top Isoliq A | Top Gyps Rid | Top Bite Rims



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

DSTVITOflexi



<b>Part no</b>	52100S	52100	52101	52102	52110S	52110	52111	52112	52120S	52120	52121	52122
<b>Color</b>	Pink	Pink	Pink	Pink	L.Pink	L.Pink	L.Pink	L.Pink	D.Pink	D.Pink	D.Pink	D.Pink
<b>Size</b>	S	M	L	XL	S	M	L	XL	S	M	L	XL
<b>Part no</b>	52140S	52140	52141	52142	52130S	52130	52131	52132	52150S	52150	52151	52152
<b>Color</b>	Pink.O	Pink.O	Pink.O	Pink.O	H.Pink	H.Pink	H.Pink	H.Pink	C.Pink	C.Pink	C.Pink	C.Pink
<b>Size</b>	S	M	L	XL	S	M	L	XL	S	M	L	XL

DSTVITOflexi is a biocompatible denture base material has exceptional aesthetic and physical properties. An excellent flexible thermoplastic is ideal for creating compliant and comfortable removable partial and mini dentures, over-dentures on implants, framework, and unilateral dentures. It resists moisture and stains. Is easy to adjust, finish and polish. Repairable and rebasable by use a fuse liquid with injection method.

- | Biocompatible
- | Essentially unbreakable
- | Aesthetic and natural look
- | High translucency
- | Comfortable and excellent Fit
- | Easy to finish and polish
- | Provides superior functionality
- | High dimensional Stability
- | Lightweight yet durable
- | Virtually invisible
- | Excellent color stability
- | Minimal Shrinkage
- | Low water absorption
- | High Stain Resistance
- | Dynamic Strength Flexibility
- | High Impact Resistance

**Package:** Box 5 pcs Ø25mm cartridge  
Bulk of 100gr, 300gr, 500gr, 1kg



DSTVITOacry



<b>Part no</b>	52200S	52200	52201	52202	52203	52210S	52210	52211	52212	52213
<b>Color</b>	Pink	Pink	Pink	Pink	Pink	L.Pink	L.Pink	L.Pink	L.Pink	L.Pink
<b>Size</b>	S	M	L	XL	XXL	S	M	L	XL	XXL
<b>Part no</b>	52230S	52230	52231	52232	52233	52220S	52220	52221	52222	52223
<b>Color</b>	Pink.O	Pink.O	Pink.O	Pink.O	Pink.O	D.Pink	D.Pink	D.Pink	D.Pink	D.Pink
<b>Size</b>	S	M	L	XL	XXL	S	M	L	XL	XXL

DSTVITOacry is a Free Monomer thermoplastic is an advanced biocompatible material for injection molding to prepare total and partial dentures, implant retained over-dentures and also suitable in combination with porcelain artificial teeth for a complete non-allergenic denture. It is the next-generation development for the replacement of conventional acryl. Due to its monomer-free composition, it does not cause allergenic reactions. Featuring exceptional properties such as high strength, lightweight, elasticity and natural translucence, it is more resistant to breakage than conventional resins.

- | Biocompatible
- | Monomer-free, hypoallergenic
- | Aesthetic, translucent and comfortable
- | Easy to finish and polish
- | Strong, durable, and lightweight
- | High standard of precision
- | Excellent retentivity
- | High impact resistant
- | High color stability
- | Minimal water absorption
- | Repairable and rebasable

**Package:** Box 5 pcs Ø25mm cartridge  
Bulk of 100gr, 300gr, 500gr, 1kg



DSTVITOlux-f



<b>Part no</b>	52000S	52000	52001	52002	52003	52010S	52010	52011	52012	52013	52020S	52020	52021	52022	52023
<b>Color</b>	Pink	Pink	Pink	Pink	Pink	L.Pink	L.Pink	L.Pink	L.Pink	L.Pink	D.Pink	D.Pink	D.Pink	D.Pink	D.Pink
<b>Size</b>	S	M	L	XL	XXL	S	M	L	XL	XXL	S	M	L	XL	XXL
<b>Part no</b>	52040S	52040	52041	52042	52043	52030S	52030	52031	52032	52033					
<b>Color</b>	Pink.O	Pink.O	Pink.O	Pink.O	Pink.O	Clear	Clear	Clear	Clear	Clear					
<b>Size</b>	S	M	L	XL	XXL	S	M	L	XL	XXL					

This innovative flexible dental resin allows users to easily create comfortable and compliant removable dentures. Created from a unique thermoplastic material, DSTVITOlux-f has exceptional aesthetic and physical properties that enable highly comfortable and natural looking partials. This rigid, yet flexible polymer is an ideal solution.

**DSTVITOlux-f** is biocompatible, high-performance polymer offers all the advantages of acrylic, enabling the creation of aesthetic, natural-looking comfortable prosthetic applications, including full and partial dentures, splints, clasps and more. Resistant to both moisture and stains, DSTVITOlux-f is easy to adjust, finish and polish. It sets faster than other materials, and because of its flexible properties, patients enjoy an improved fit.

- | Revolutionary, Rigid yet flexible denture base material
- | Easily repaired and relined by conventional acrylic
- | luxury comfort, durability and moderate flexibility
- | Strong yet flexible polymer is easy to work with
- | Repels plaque and odor
- | Exceptional aesthetic and physical properties
- | Enables a more perfect fit

**Package:** Box 5 pcs Ø25mm cartridge  
Bulk of 100gr, 300gr, 500gr, 1kg





## DSTVITOkeep



<b>Part no</b>	52400	52401	52440	52441
<b>Color</b>	Tan	Tan	Pink	Pink
<b>Amount</b>	50gr	100gr	50gr	100gr

DSTVITOkeep is a semi-crystalline thermoplastic with excellent mechanical and chemical-resistance properties. As a high performance polymer is the best alternative for non-metallic restorations to fabricate full and partial removable dentures and implant overdentures frameworks, crown and bridge construction.

- Novel biomaterial as an alternative to zirconia
- Individual adjustment and prevents chipping
- Optimal work flow
- Allows a high standard of precision and finishing
- Perfect aesthetic results
- Magnificent feeling
- High resistance to abrasion
- Ideal mechanical properties
- Particularly attractive
- Biocompatibility, Greater Versatility





Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## DSTVITOp last

Softens in hot water



<b>Part no</b>	52700	52701	52702	52710	52711	52712
<b>Color</b>	Blue	Blue	Blue	Yellow	Yellow	Yellow
<b>Amount</b>	300 gr	500 gr	1 kg	300 gr	500 gr	1 kg
<b>Part no</b>	52720	52721	52722	52730	52731	52732
<b>Color</b>	Pink	Pink	Pink	Green	Green	Green
<b>Amount</b>	300 gr	500 gr	1 kg	300 gr	500 gr	1 kg

DSTVITOp last is an innovative alternative for conventional custom tray material. Designed for a variety bite registration -fast and easy to use. Friendly and convenient for the patient and dentist - comfort custom trays, bite registration, fore and side bite detection in 1.2.3. Features good working time along with excellent malleability and adaptation to model. Easy trimming characteristics reduce effort.

- Affordable, cost effective
- Biocompatible and non-allergenic
- Newest alternative for conventional tray material
- Softens in hot water
- Rigid yet unbreakable
- Comfortably handled and shape
- Can be simply cut and fit in softened state
- Extraordinary accuracy
- No odor, very durable



DSTVITOCetal



<b>Part no</b>	52300S	52300	52301	52302	52310S	52310	52311	52312	52320S	52320	52321	52322
<b>Color</b>	Pink	Pink	Pink	Pink	A2	A2	A2	A2	A3	A3	A3	A3
<b>Amount</b>	S	M	L	XL	S	M	L	XL	S	M	L	XL
<b>Part no</b>	52330S	52330	52331	52332	52340S	52340	52341	52342				
<b>Color</b>	A3.5	A3.5	A3.5	A3.5	C-D	C-D	C-D	C-D				
<b>Amount</b>	S	M	L	XL	S	M	L	XL				

DSTVITOCetal is a innovative biocompatible, non-allergenic flexible thermoplastic resin - It's ideal metal replacement aesthetic solution for partial denture frameworks, aesthetic clasps, crowns, bridges and other metal-free removable applications by injection molding systems in combination with our DSTVITOLine of thermoplastics and conventional acrylics.

DSTVITOCetal popularity due to excellent physical and mechanical properties and has fostered use for creating aesthetic dental appliances with flexibility, stability and comfort.

Low modulus of elasticity, which allow for its use in larger retentive undercuts this, may be advantageous in clinical situations in which aesthetics or periodontal health are priorities.

- | Aesthetic and comfortable
- | Very low weight
- | High tensile and flexural strength
- | Superb dimensional stability
- | High elastic memory
- | Fatigue resistance
- | Optimal rigidity
- | Excellent impact resistance
- | Easy to polish and easy to clean

**Package:** Box 5 pcs Ø25mm cartridge  
Bulk of 100gr, 300gr, 500gr, 1kg



# DST –FLEXI DENTURE MATERIALS

## DST CLEAN & RENEW

The DST Clean/Renew concentrated denture cleaner powder, is an excellent solution and recommended for cleaning thermoplastic and acrylic dentures.



- | Do not apply while denture is in mouth, use as instructed with dentures removed.
- | Can be mixed and used for up to 7 days before being replaced with a new mix.
- | One spoonful of powder is mixed with 7 oz. of warm water.

Part num	Color	Amount
52-690	powder	4 oz

## DST LUSTERLIQ

DST LUSTERLIQ can be used to the finish gloss on any acrylic and thermoplastic denture appliances.



Part num	Color	Amount
52-620/4	Transparent	4 oz
52-620/8	Transparent	8 oz

- | Easy to apply
- | resulting in the quickly, high-gloss luster.
- | Provides a brilliant surface finish
- | Leaves a pleasant fresh taste

## DST LUBRILIQ F

DST Lubriliq it is perfect for lubricating aluminum cartridges prior to insertion into a press oven, or for greasing the flasks.



Part num	Color	Amount
52-610/4	Transparent	4 oz
52-610/6	Transparent	6 oz
52-610/8	Transparent	8 oz

- | Multi-purpose release lubricant
- | easy and safe release
- | Suitable for all injection molding machines
- | Easy to application assures even coverage



Dental Materials Production

www.dst-dental.com | Email: info@dst-dental.com

## DST ISOLIQ F



**REVOLUTIONARY NEW FORMULA**  
**Deflask and get smooth, glaze shine effect**

DST Isoliq F is a extra-premium separator specially formulated for thermoplastic denture injections molding. A single, thin coat will provide a smooth and special glaze shine effect after injection. Best results - injections right out of the flask, saving time and money in finishing and polishing.

- | To get a perfect results
- | Smooth, glaze shine effect
- | Thin isolation layer
- | Fast drying
- | Water repellent
- | High heat resistance

Part num	Color	Amount
52-600/4	Transparent	4 oz
52-601/6	Transparent	6 oz

## DST CEM F



Our DST Cem F is a very strong transparent glue that bonds artificial acrylic teeth into the stone or the flask securely and firmly, after boil out - prior to injection. DST Cem F is mild enough to keep from bonding the teeth to the stone post-injection.

- | Easy use with applicator of thin liquid glue.
- | Fast drying without residue.
- | Will not stain or damage surface of teeth.

Part num	Amount
52-650	20 ml


## DST FUSE FLEX



This is a liquid fusing bond that is formulated for repairs, linings or additions to polyimide resin appliances and partials. DST Fuse Flex begins to deteriorate existing bonds, and creates a new bond when injected into the area. The new bond that is created is stronger and virtually undetectable.

- | Makes repairs and additions to existing partials immediately prior to re-injection.
- | Use for minor repairs and additions to existing DST Unbreakable Partial
- | Use only when flasking and injection are absolutely necessary.
- | Not to be used in the mouth. Must be used in dental laboratory

Part num	Amount
52-640	1 oz

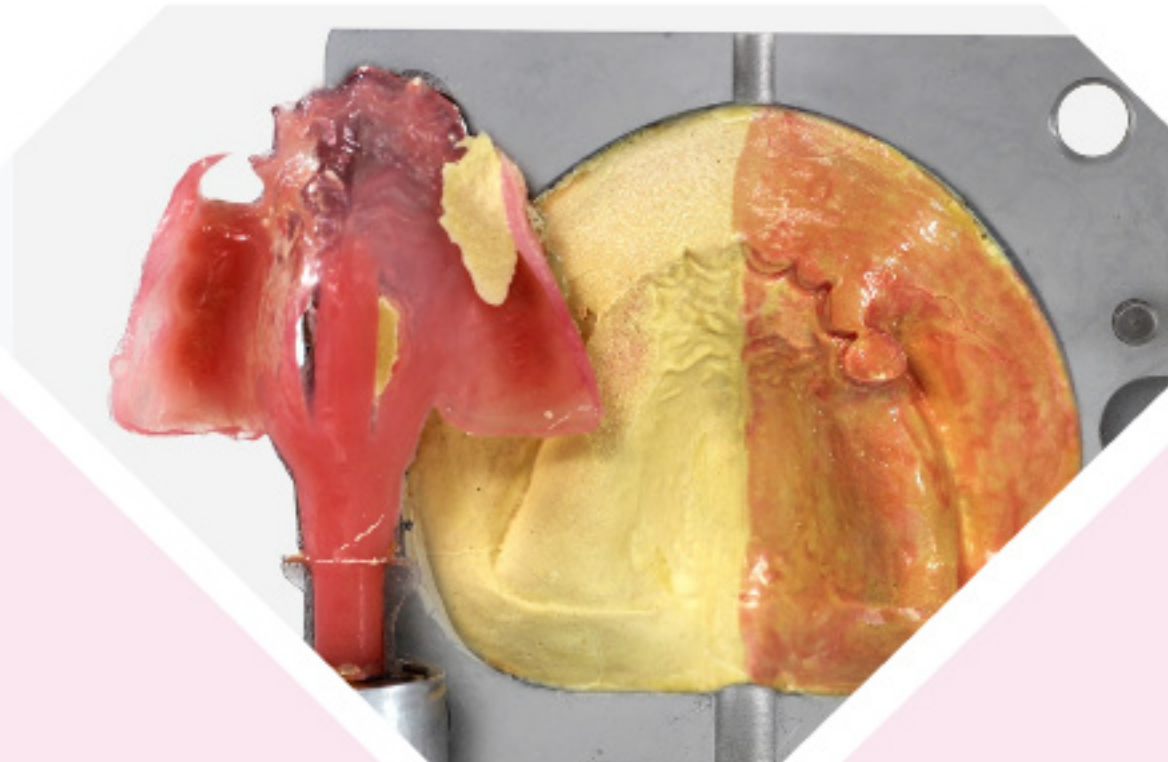

## Top Isiliq A



### PLASTER & ACRYLIC SEPARATOR


Tin foil substitute for separation of acrylic appliances from stone models

TOP ISOLIQA A is a alginate-based separating solution formulated for use between gypsum stone & plaster and acrylic resin . For the insulation of plaster surfaces against denture resins (heat and cold curing) and composites.

- | Prevents cured acrylic denture from sticking to stone plaster cast.
- | Thin separating film.
- | Short action time and no formation of layers
- | Achieving smooth surfaces
- | Optimum separation - High quality
- | Works with cold or hot models

Part num	Color	Amount
26-750/16	Pink	16oz
26-750/32	Pink	32oz

## Top Gyps-Rid

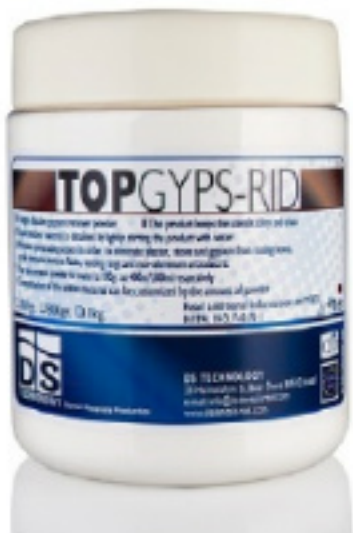

- | Highly effective gypsum remover in powder form.
- | Becomes fully efficient in a 35-45 C° ultrasonic bath.
- | Stirring the powder with water creates a solution to the viscosity required.
- | Dilution ratio is 100g powder to 400-500ml water.
- | Keeps the utensils shiny and clean.
- | Active material concentration can be customized by powder/water proportions.

Part num	Color	Amount
26-500	Powder	500 g

## Top Bite Rims

Anatomically formed, flexible and comfortable, our Top Bite Rims Wax is able to achieve exact bite registration in one application. These pre-formed rims include a basal channel for good fit and registration, and come in a hardness that accepts excellent tooth detail, yet is comfortable to bite.



Part num	Color	Amount
26-130	Red-flavored preformed	50 pcs

- | Flavored and comfortable for patient.
- | Standard jaw sizes.
- | Contains 50 pieces per box.









# Ceramic Products

Top Ceraliq | Top Cera Glazer | Top Marker | Top Cerafix Gel  
Top Cera Contour | Top Cera Sep | Top Cera Sep - Pen Applicator  
Top Cera Sep Spray



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## Top Ceraliq



Our Top Ceraliq is a premier condensing liquid that assists any brand of porcelain powder to form a firmer, easier to control, hard finish. This product enhances the hardening properties, creating an ionic charge to automatically condense individual porcelain particles, helping them become denser and a tighter unit, enhancing the color, longevity and predictability of the application while eliminating air bubbles and particle clumping.

- | Longer moisture phase for predictable, reliable results.
- | No fear of dry-out, no need for frequent re-wetting and dilution.
- | Fast build-up, minimal shrinkage and more efficient condensing.
- | Shorter working time and a stronger, more natural looking result.

Part num	Color	Amount
26-510/4	Transparent	4 oz
26-510/16	Transparent	16 oz
26-510/32	Transparent	32 oz

## Top Cera Glazer



Our Top Cera Glaze is a specially formulated liquid glaze medium that, when utilized during glaze firing results in a natural looking cosmetic restoration.

- | Universal Top Cera Glaze can be easily mixed with Hi Cera Solution for exact match.
- | Mix Top Cera Glaze liquid with glaze or stain powders for smooth consistency.

Part num	Color	Amount
26-520	Transparent	15 ml
26-521	Transparent	30 ml
26-520/4	Transparent	4 oz

## Top Marker



A fast-drying low-viscosity marking pigment for detecting high spots and tight spots when seating castings, Top Marker provides high-contrast, accurate smears transferred at the slightest pressure on the opposing side to identify 'problem areas' on interproximal and occlusal areas.

- | Applies with brush.
- | Will not stain or discolor, does not contaminate metal.
- | Easy removal with water and soft brush.

Part num	Color	Amount
26-940	Blue	1/2 oz
26-941	Yellow	1/2 oz
<b>Thinner</b>		
26-945/1	Transparent	1/2 oz



## Top Cerafix Gel

High quality holding gel for porcelain veneers & full-ceramic restorations on the die  
 Top Cerafix Gel is a high quality temporary adhesive gel for holding crowns, Zirconia copings, veneers, wax and small metal castings in place during fabrication - procedures the build-up and contour processes.  
 Carefully formulated with excellent hold veneers power, this gel is the perfect 'third hand' helper for temporary fitting.

- | Easy to Use
- | Top quality gel
- | Is especially useful in holding Pressed Veneers
- | Excellent holding and easy separation
- | Easily remove restoration from die
- | Burns without residue
- | Water-soluble
- | Transparent
- | Saves Time and Money

Part num	Color	Amount
26-515	Transparent	4 ml
26-516	Transparent	10 ml



## Top Cera Contour



Top Cera Contour is a surface structures indicator. Ready for use by brush. Easily applied a thin layer. Visualization of imperfections, flaws, defects, irregularities allowing for the perfect match. Now you can achieve amazing results restorations in a short time and easily.

- | Defects and imperfection becomes visible
- | Simply and neatly smears
- | Dries within seconds
- | Easy clean with steam
- | Providing you the ability to create perfect restorations

Part num	Color	Amount
26-526	Gold	1/2 oz
26-527	Silver	1/2 oz



## Top Cera Sep

Our Top Cera Sep is perfect for sealing a die or model against the absorption of moisture when molding porcelain. Enhances release properties and prevents sticking of porcelain when separating from the die, prevents delicate, intricate parts from remaining with the die. Top Cera Sep does not cause discoloration of porcelain product, and burns out clean.



- | Separates porcelain from model material on gingival surface of pontics, in marginal porcelain regions and in adjacent contact zones.
- | Thin layers prevent end product faults.

Part num	Color	Amount
26-530	Transparent	1oz

### Top Cera Sep Pen Applicator

Pen Applicator- For the separation of plaster from porcelain and wax  
Recommended use on the sealed or hardened surface area

- | Reliable separating efficiency.
- | Easy application due to pen shape.
- | Easy and quick removal of the porcelain set-up
- | Prevents the moist porcelain to stick to plaster.
- | Enables a uniform and thin application.
- | Free of solvents.

Part num	Color	Amount
26-531	Transparent	1/2 oz

Refill

26-532/4	Transparent	4 oz
26-532/8	Transparent	8 oz



## Top Cera Sep Spray



Porcelain Release/Separating Agent.  
Prevents porcelain from sticking to the stone model at the build-up.

- | Enhances release and prevents sticking of porcelain.
- | For release of onlays, occlusal & proximal contacts & resealing margins.
- | Conveniently packaged product - Pump Spray.
- | Easy to apply.
- | Does not cause discoloration.
- | Burns out clean.

Part num	Color	Amount
26-535/4	Transparent	4 oz
26-535/8	Transparent	8 oz



Dental Materials Production

www.dst-dental.com | Email: info@dst-dental.com



# Model Production

TOPDEB Debubblizer | Top Vacudeb | Top Sep | Top Smart Sep  
 Top shield&Shine Spray | Top Blockout Wax | Top Die Hardener& Sealer  
 Top Die Hardener& Sealer - Pen Applicator | Top Die Sealer



Dental Materials Production  
[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## TOPDEB Debubbler



When pouring into molds for casting dental models, surface tension can push back the casting material and form imperfections and air bubbles. Our Top-Deb Debubbler prevents this by stabilizing the surface of the mold, allowing material to flow into the thinnest areas and adhere completely to the mold, allowing accurate, highly detailed and dense castings. It can also be used in treating stored impressions and wax patterns, to aid in removal of dust particles and oily residues which accumulate over time.

Universal application, including silicone, rubber base or alginate impression material. Transparent, blue tint. Includes refillable misting spray container for easy application.

Part num	Color	Amount
26-400/8	Blue - Transparent	8 oz
26-400/16	Blue - Transparent	16 oz
26-400/32	Blue - Transparent	32 oz

## Top Vacudeb



A surface tension reducer, Top Vacudeb allows an attachment to achieve close adaption and maximum contact with all surfaces of a wax pattern and with all impressions prior to investing, without trapping or forming air bubbles or voids on the surface of models and castings.

By managing surface tension, an improved flow rate of investment material and strict adherence to the impression ensures that the casting is reproduced perfectly and accurately.

- May be used under vacuum.
- Fast drying increases productivity.
- Eliminates air entrapments, bubbles and voids.
- Facilitates fluid flow to produce sharper detail.
- Supports silicone, wax, rubber and alginate patterns.
- Transparent pink tone.

Part num	Color	Amount
26-220/8	Pink- Transparent	8 oz
26-220/16	Pink- Transparent	16 oz
26-220/32	Pink- Transparent	32 oz

## Top Sep



A non-oily anti-sticking separating agent, and dries on contact. Since it is alcohol-based, Top Sep quickly dries to a hardened finish on the treated surface, allowing immediate pouring of investment into dies or molds. There is no film left behind.

- Environmentally safe; contains no silicones or chlorofluorocarbons.
- Soaks into stone or plaster to seal pores and assure smooth finish of investment.
- Hardens and smooths surface of stone or plaster, and makes it water-proof.
- Quick cleanup if sprayed on trays, vessels, flasks, sinks, etc.
- Leaves no film.

Part num	Color	Amount
26-200/8	Yellow Transparent	8 oz
26-200/16	Yellow Transparent	16 oz
26-200/32	Yellow Transparent	32 oz

## Top Smart Sep

TOP SMART SEP specially formulated for separating casting and modeling waxes from porcelain, stone, composite, cement, acrylic resins and metal.



- | Universal application.
- | Easy to apply.
- | Restorations exhibit smoother and cleaner surfaces.

Part num	Color	Amount
26-201/8	Transparent	8oz
26-201/16	Transparent	16oz
26-201/32	Transparent	32oz

## Top Shield & Shine

High gloss coating to all gypsum models - from orthodontic to the practical models



- | Fast drying
- | Easy to apply
- | Get the professional shine luster and hardens
- | Prevents staining & chipping of models
- | Provides a smooth, glossy stone models
- | Ideal from orthodontic models
- | Professional quality look when presenting models - improves case presentations.
- | Works on all dental stones & plasters
- | Thickness of less than one micron
- | Accuracy & quality are maintained
- | Not create occlusal interference.

Part num	Color	Amount
26-925/4	Transparent	4oz
26-925/8	Transparent	8oz
26-925/16	Transparent	16oz
26-925/32	Transparent	32oz

## Top Blockout Wax

A very hard, high temperature wax, our Top Blockout Wax is designed to fill undercuts and process induced cavities, cover sharp edges and smooth out rough surfaces of dies prior to molding attachments. Excellent carvability and very good scraping properties.



- | Fast and secure application to undercuts on dies for molding of comfortable attachments.
- | Adapts to die contours while soothing.
- | Does not burn without residue.
- | Will not coalesce with dipping wax, cervical wax or other modeling waxes due to very high melting point.
- | Does not discolor plaster model after boil-out.
- | For blocking out bite guard splints and dies.
- | 105 °C melting point.

Part num	Color	Amount
26-500	White Transparent	42/45g
26-501	Pink Transparent	42/45g



## Die Hardener / Sealer

Our Die Hardener is an absorbing resin that seals, hardens and strengthens stone models and dies. As it is absorbed into the stone surface, it permeates into all of the fissures and microscopic pores, bonding the entire stone by fusing the particles together. It resists chipping by leaving a zero micron infusion to the surface layer bonded to the resin permeating the stone, assuring an accurate and durable mold.



- | 2 in 1
- | Seal pores and high bond strength
- | Fast drying to all dental stones
- | Resistant against separating agent
- | Waterproof, abrasion-resistant surface
- | Protects from chipping and cracking
- | Produces a smooth, dense, hard surface
- | Seal And hardness dies and stone models

Part num	Color	Amount
26-935	Transparent	1 oz
26-935/4	Transparent	4 oz
26-935/8	Transparent	8 oz

Thinner

26-911	Transparent	1 oz
--------	-------------	------

Special thinner for setting preferred consistency Especially adapted to the texture of the TOP DIE HARDENER & SEALER

### Die Hardener / Sealer Pen Applicator

- | Reliable separating efficiency.
- | Comfortable and Easy application due to pen shape.
- | Economical and on-the-spot application;
- | Fast drying in only seconds. Hardens, seals, smoothes
- | Fills open micro pores on the plaster surface
- | Abrasion-resistant surface. Dirt&water-repellent
- | Resists chipping and cracking.

Part num	Color	Amount
26-931	Transparent	1/2 oz

Refill

26-931/4	Transparent	4 oz
26-931/8	Transparent	8 oz



## Top Die Sealer

Our Top Die Sealer is brush applied to easily seal and water-proof stone die and model surfaces. It protects margins and reduces abrasion during operations, protecting the die from chipping, breaking or sticking. It also can provide an excellent base layer for our Top Cera Sep separating substance between porcelain and cast. It separates porcelain from model material on gingival surface of pontics, in both marginal porcelain regions as well as adjacent contact zones.



- | Protects preparation borders of crowns and inlays.
- | Prepares porcelain shoulder layers.
- | Adds negligible film thickness.
- | Thin consistency allows penetration into stone for firm seating.
- | Transparent surface becomes hard and water-proof.
- | Requires only one application, although additional thin layers can be applied after first coat is dry.

Part num	Color	Amount
26-920	Transparent	1 oz

Thinner

26-911	Transparent	1 oz
--------	-------------	------







# Die Preparation

Top PC Spacer | Top PC Spacer - PCS ISO | Top PC Spacer - KITS  
Top Die Spacer | Top Die Lube | Top Die Lube - Pen Applicator  
Top Rubbsep | Top Blockout Wax



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## Top PC Spacer

### FOR PRESSABLE CERAMIC TECHNIQUE

For the creation of highly aesthetic all-ceramic restorations, including inlays and onlays, veneers and partial crowns. Enhances the visual perfection by emulating the natural shade. Its natural tint provides optimum compatibility with an accurate shade matching the patient's dentition to perfect results

**Applies quickly and evenly , Dries quickly , Smooth surface, Simulating accurate die cement spacing.**

**ATTENTION:** Use only with PCS ISO - special separating agent

The use of incompatible isolating agents could provoke unforeseen reactions

**Available in :** 9 IPS® Natural Die Material® shades \* and milk shade, 1/2oz

\*IPS Natural Die Material® is a registered trademark of Ivoclar Vivadent.



<b>Part no</b>	52-2801	52-2802	52-2803	52-2804	52-2805	52-2806	52-2807	52-2808	52-2809	52-2810
<b>Color</b>	PCS 1	PCS 2	PCS 3	PCS 4	PCS 5	PCS 6	PCS 7	PCS 8	PCS 9	PCS Milk
<b>Amount</b>	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz	1/2 oz

**Top PCS Thinner - Part no -52-285, Color - Transparent, Amount - 1/2 oz**

## PCS ISO

PCS ISO is specially formulated as a compatible preparation with Top PC Spacer and recommended to be used together. The use of incompatible isolating agents could provoke unforeseen reactions. Forms a very thin layer which prevents surface falsification. PCS ISO is formulated from environmentally safe materials, contains no solvents and will not attack die varnishes or waxes. It won't evaporate either.

ADVANTAGE IS a optimum compatibility and NOT damage the - TOP PC SPACER- TOOTH COLOURED DIE SPACERS



- | Facilitates the removal of wax patterns .
- | Prevents the die material from sticking to the ceramic material.
- | Ensures precise and comfortable fit of wax couplings.
- | Developed for use with our Top PC Spacer for speedier and simple application.
- | Low surface tension

Part num	Color	Amount
52-284	transparent	1/2 oz

Top PC Spacer -KITS



COMPARISON CHART

Top PC Spacer Shades	PCS1	PCS2	PCS3	PCS4	PCS5	PCS6	PCS7	PCS8	PCS9	PCS MILK
IPS® Die Natural Material Shades	ND1	ND2	ND3	ND4	ND5	ND6	ND7	ND8	ND9	
Vita®	A1,A2	A1,A2 A3, B2	B3,B4	A3.5	A4		C1,C2, C3, D1, D2, D3			

COMPLETE KIT

Part num - 52-280  
 All 10 Natural Shades in 1/2 oz bottle: PCS 1, PCS 2, PCS 3, PCS 4, PCS 5, PCS 6, PCS 7, PCS8, PCS 9, PCS MILK  
 PCS Thinner 1/2 oz  
 PCS ISO 1/2 oz

STARTER KIT

Part num - 52-281  
 6 Natural Shades in 1/2 oz bottle: PCS 2, PCS 3, PCS 5, PCS 7, PCS 9, PCS MILK  
 PCS Thinner 1/2 oz  
 PCS ISO 1/2 oz



Dental Materials Production

www.dst-dental.com | Email: info@dst-dental.com

## Top Die Spacer



In order to control the relief and consistency of dies and castings, a uniform spacing before wax-up of a mold is essential. Our Die Spacer is a precise paint-on liquid that maintains a 12 micron homogeneous thickness, flows on smoothly and bonds to the gypsum or stone mold without running or clumping.

- | Precise control and viscosity to allow uniform coatings.
- | Fast drying, allowing quick application of multiple coats.
- | Contrasting color for definition and inspection.
- | Outstanding resistance to pressurized steam, chipping or process damage.

Part num	Color	Amount
26-900	Blue	1 oz
26-900	Blue	1/2 oz
26-901	Red	1/2 oz
26-902	Gold	1/2 oz
26-903	Silver	1/2 oz
Thinner		
26-910	Transparent	1 oz

## Top Die Lube

Top Die Lube is a fast-drying, non-oily, residue-free release film that is easy to apply to a die, providing a wet surface liner for wax castings. The thin finish allows for greater accuracy of the casting, no pattern distortion and clean separation, even after standing overnight. Excellent characteristics and performance off of stone, plated or epoxy dies.



- | Can be used with dripping wax or regular wax.
- | Easy removal, even after time, with no residue attached to casting or die.
- | Die spacer will not be removed with casting.
- | Contains alcohol, green color.

Part num	Color	Amount
26-800	Green	1 oz

### Top Die Lube Pen Applicator

Pen Applicator -For the separation of plaster from wax

- | Reliable separating efficiency.
- | especially for the wax dipping technique.
- | Easy and quick removal of the wax-up from the model without distortions.
- | Economical Easy application due to pen shape.
- | Enables a uniform and thin application.
- | Adjusted to our Top waxes and Spacers.
- | Can be used with dripping wax or regular wax.



Part num	Color	Amount
26-801		1/2 oz
REFILL		
26-802	Transparent	100ml

## Die Preparation

### Top Rubbsep



to apply rubber film that forms a peel-able spacer between the die and the attachment during

- Provides accurate spacing for precision attachments.
- Rubber coating prevents free monomer of resin from curing teeth.
- Provides thermal protection for teeth.
- Prepares impression pour release agent for crowns.

Part num	Color	Amount
26-960	Pink	1/2 oz
26-961	White	1/2 oz



- Fast and secure application to undercuts on dies for molding of comfortable attachments.
- Adapts to die contours while soothing.
- Does not burn without residue.
- Will not coalesce with dipping wax, cervical wax or other modeling waxes due to very high melting point.
- Does not discolor plaster model after boil-out.
- For blocking out bite guard splints and dies.
- 105 °C melting point.

Part num	Color	Amount
26-500	White Transparent	42/45g
26-501	Pink Transparent	42/45g



Dental Materials Production

www.dst-dental.com | Email: info@dst-dental.com



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)



# Cad/Cam

Top Scan Wax | Bari x-ray



Dental Materials Production  
[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## Top Scan Wax

### CAD/CAM opaque

Our Top Scan Wax for CAD/CAM is a specially formulated laser-opaque and light-opaque wax for sculpturing according to the CAD/CAM procedure. Guarantees superb scanning results which enhances exact fitment. It builds up fast, sets rapidly and is quickly available for carving and scraping operations.



- | Highly opaque for prime scan operation.
- | No powder spray necessary.
- | High stability after hardening maintains dimensional integrity.
- | Superb modeling qualities.
- Not to be used in casting.

Part num	Color	Amount
26-550	Beige- opaque	42-45 g
26-551	purple - opaque	42-45 g
26-552	Turquoise- opaque	42-45 g

## Bari x-ray

To mix with acrylics for X-ray feature, for optimum diagnostic imaging.

Bari x-ray is excellent for adding and mixing with Acrylic Resin, provides a high contrast identifier in X-Ray scanning pictures with consistent texture, for the fabrication of methacrylate - based acrylic appliances and used wherever a radiopaque acrylic dental device is useful in establishing an anatomical location or site , such as a provisional implant stent.



- | Ideal for production of radiopaque acrylic dental device
- | diagnostics, and planning of treatments.
- | Identifies primary implant location.
- | borders are sharp for precise implant planning

Part num	Amount
52-700	4 oz





# Wax Products

Dipping Waxes | Crown and Bridge Waxes | Special Waxes  
Sticky Waxes | Top Wax Finish | Uni-Wax Adhesive



Dental Materials Production

[www.dst-dental.com](http://www.dst-dental.com) | Email: [info@dst-dental.com](mailto:info@dst-dental.com)

## Dipping Waxes

Top Dip Wax is used in electrical dripping pots. Perfect measures allow precise and consistent results. Once liquefied, provides excellent elasticity, adaption and detail. Retains strong edges with a memory enhanced elastic feel, clean burn out and minimum shrinkage for consistent dimensional integrity.

- | Dries quickly.
- | Eliminates resealing.
- | Maximum hygroscopic expansion for casting fit requiring no additional adjustment.
- | Contains natural resins for stability and elasticity.
- | Excellent scraping and cutting properties allow best fittings and easy crown adaption.
- | Resists warping, deforming or breaking of copings.
- | Can be held in molten state for long periods of time.
- | Easy clean up.

### Top Dip Wax S



Part num	Color	Amount
26-100D	Red	100g

- | Homogeneous wax form model results.
- | Leaves no wax residue on model.
- | Retains strong edges.

### Top Dip Wax M



Part num	Color	Amount
26-200D	Orange	100g
26-210D	Pink	100g

- | Perfect elasticity properties.
- | Recommended hardness for all applications.

### Top Dip Wax H



Part num	Color	Amount
26-300D	Green	100g

- | Great for making secondary forms with conus or telescope technique.
- | Precise and consistent results.

## Crown and Bridge Waxes

### Top Modeling Wax

Our Top Modeling Wax features optimum properties for carving and scraping, due to non-elasticity, hardness without brittleness, and high surface tension. This wax has quick application and cooling properties, speeding up the work flow when building bridges, crowns and inlays.



- | Optimal raw materials allow additive wax-up techniques.
- | High density molecular structure is easily carvable and leaves smooth surfaces.
- | Strong enough for elaborating detailed cusp tips.
- | Low shrinkage, low expansion, opaque.
- | Clean burnout, ash free, leaves no traces.

<b>Part no</b>	26-200	26-201	26-202	26-203	26-204	26-205	26-206
<b>Color</b>	Gray- opaque	Blue- opaque	Green- opaque	Violet- opaque	Brown- opaque	Red- opaque	Beige - opaque
<b>Amount</b>	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g

### Top PC Modeling Wax

Specifically formulated for press-able ceramic systems, our Top PC Modeling Wax is filtered several times to completely eliminate any opaque particles which might negatively react with the investment material or casted object. Like our Top Modeling Wax, this PC version retains excellent sculpting, caving and scraping characteristics due to our exclusive high density molecular structure.



- | No color contamination.
- | Ash free for clean burnout.
- | Compatible with all press-able systems.
- | Low shrinkage.
- | Fast application and quick cooling.
- | High surface tension for precise carving and smooth finish.
- | Ash Free

<b>Part no</b>	26-400	26-401	26-402	26-403	26-404	26-405	26-406
<b>Color</b>	Gray	Blue	Green	Violet	Yellow	Bordeaux	Beige
<b>Amount</b>	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g

### Top Sculpturing Wax

Offering a very firm consistency that allows excellent handling and carving suitability, our Top Sculpturing Wax will not become brittle, crack or fracture during process activities or curing. This is a harder wax that is easy to scrape, carve and create detailed occlusal anatomy. When heated and in pouring or dipping form, the viscosity is easily controlled, covering large areas expediently and efficiently.



- | Low contraction factor allows precision fit.
- | Low shrinkage for unsurpassed fabrication accuracy.
- | Clean burnout due to resin and ash free manufacture.
- | Stable inner elasticity allows for true reproduction of occlusal stops and contact points.
- | Ash Free

<b>Part no</b>	26-100	26-101	26-102	26-103	26-104	26-105	26-106
<b>Color</b>	Gray	Blue	Neon Green	Neon Yellow	Orange	White	Beige
<b>Amount</b>	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g	42-45 g

Special Waxes

Top Diagnostic Wax

Specially formulated for diagnostic use, our Top Diagnostic Wax has an elevated melting point and a natural finish for lifelike results. Modeling properties match the best properties of our leading sculptural waxes while providing extra elasticity for accurate representation of detail, surface and finish.



Part num	Color	Amount
26-450	Dentine-opaque	42-45 g
26-451	Write-transparent	42-45 g
26-452	Write- opaque	42-45 g

- | Extremely stable, low in contraction.
- | Ideal for elaborating wax-ups.
- | Easy to scrape and carve.
- | Color corresponds to tooth enamel.
- | Not to be used in casting.

Top Cervical Wax

Ensuring a perfect marginal seal for inlays, onlays and crowns requires our completely shrink-free Top Cervical Wax. Adheres perfectly and comfortably to the coping wax, its smooth flow provides precise margins and secure fitment. Carving and scraping touch up can be carried out without tearing, allowing construction of very thin margins.



Part num	Color	Amount
26-350	Red	42-45 g

- | Non-rebounding and tension-free.
- | No pulling or flaring.
- | Compatible with Top Dip Wax and Top Sculpturing Wax.
- | Easy flowing, non-smear qualities allow precision application and a sharp and accurate margin.
- | Ash free.



Dental Materials Production

www.dst-dental.com | Email: info@dst-dental.com

## Special Waxes

### Top Scan Wax

Our Top Scan Wax for CAD/CAM is a specially formulated laser-opaque and light-opaque wax for sculpturing according to the CAD/CAM procedure. Guarantees superb scanning results which enhances exact fitment. It builds up fast, sets rapidly and is quickly available for carving and scraping operations.



- | Highly opaque for prime scan operation.
- | No powder spray necessary.
- | High stability after hardening maintains dimensional integrity.
- | Superb modeling qualities.
- | Not to be used in casting.

Part num	Color	Amount
26-550	Beige- opaque	42-45 g
26-551	purple - opaque	42-45 g
26-552	Turquoise- opaque	42-45 g

### Top Blockout Wax

A very hard, high temperature wax, our Top Blockout Wax is designed to fill undercuts and process induced cavities, cover sharp edges and smooth out rough surfaces of dies prior to molding attachments. Excellent carvability and very good scraping properties.



- | Fast and secure application to undercuts on dies for molding of comfortable attachments.
- | Adapts to die contours while soothing.
- | Does not burn without residue.
- | Will not coalesce with dipping wax, cervical wax or other modeling waxes due to very high melting point.
- | Does not discolor plaster model after boil-out.
- | For blocking out bite guard splints and dies.
- | 105 °C melting point.

Part num	Color	Amount
26-511	Light Pink- opaque	42-45 g
26-501	Pink Transparent	42-45 g
26-500	white opaque	42-45 g
26-510	white Transparent	42-45 g

### Top Wax Finish

Our Wax Finish is a premium smoothing wax sealer, applied to the occlusal surface of a modellation; this product smooths out potential irregularities and levels small scratches and unevenness. Floods any fissure depths that can be inaccessible for a fissure frill.



- | Creates a very smooth, monogeneous surface.
- | Dries in 2-3 minutes, depending on thickness of layer.
- | Finer surface after modelling reduces elaboration time by up to 50%.
- | Occlusal and approximal contact points remain constant after polishing.
- | Reinforcement of contact points.

Part num	Color	Amount
26-000	Transparent	1/2 oz
Thinner		
26-005	Transparent	1/2 oz

### Uni-Wax Adhesive

For use with preformed wax patterns, our Uni Wax Adhesive is a residue free bonding agent for use with preformed wax patterns used with unhardened investment duplicate models. Adhesive can be applied thinly to ensure retention of patterns on the duplicate model or the investment model without any gap or airspace, providing a secure adherence. Thin application allows part of the liquid to soak into the investment to strengthen the bond.



- | Useful for wax and resin components.
- | Contains no heavy metals.
- | Leaves no film or residue.
- | Burns without residue.

Part num	Color	Amount
26-050	Transparent	1/2 oz

## Sticky Waxes

### Top Sticky Wax

Adhesive power is one of the remarkable properties of our Top Sticky Wax -making it the best choice for break repairs, fixing teeth, front walls and solder molds. Quick setting and high strength prevent distortion and slippage during work operations. It is essentially a 'wax glue', very tacky with low flexibility after setting, making it a top choice for positioning, holding in place and making sprue connections.



- | 'Wax glue' properties for positioning.
- | Stays viscous to make attachments.
- | High stability after hardening maintains dimensional integrity.
- | Can be cleaned up or removed with steam jet.
- | Leaves no residue or contamination.

Part num	Color	Amount
26-250	Neon-Pink	42-45g

### Top Spec Sticky Wax

Top Spec Sticky Wax with extremely strong adhesive properties, is formulated to meet the needs of professionals. It quickly and accurately fixates acrylic teeth within a silicone mould. This product was explicitly developed by our specialists and allows the technician to fixate acrylic teeth without the use of glues, assuring perfect results.



- | Fixates acrylic teeth - no adhesive application.
- | Prevents damage to the anatomical shape of acrylic teeth.
- | Quick, easy and safe to apply.
- | Extreme hardness, distortion is impossible.
- | Easily removed with steamers.
- | Leaves no residue.

Part num	Color	Amount
26-300	Red	42-45g





“Love of beauty is taste.  
The creation of beauty is art.” Ralph Waldo Emerson





## DS TECHNOLOGY LTD

29 Omarim st., Omer Industrial Park, 8496500, Israel  
info@dsdentalisrael.com | www.dsdentalisrael.com  
Tel. +972-77-2006636 Fax. +972-72-2829896

